GARANTIA / GUARANTEE

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All products manufactured by Felisa are guaranteed for one year from date of shipment from Felisa plant. All those products returned within one year will be build or replaced under the guarantee regardless of reason for failure. Electric motors, are guaranteed according to the manufacturer policies.

Exceptions may be made by Felisa on particular applications, where experience has indicated conditions are so unusual that premature failure can be expected. Transportation charges, in all cases will be at customer expense. Maximum liability is in no case to exceed the value of the Felisa product involved.

Felisa has the right to change or modify the products in order to improve their presentation and/or operation.

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The care you take in reading and following this instructions will probably determinate the satisfactory service you will receive from your pump.

NOTE: This equipment must be used only for its intended application, any alteration or modification will void the guarantee.

UNPACKING.

Carefully remove the pump from the shipping case. Preserve all paper work for future reference. If damage has occurred from shipment a claim must be filed with the carrier immediately; preserve the shipping container for inspection by the carrier. Contact your dealer or Elisa.

INSTALATION.

All Felisa vacuum pumps should be mounted in a horizontal plane. Rubber bumpers are supplied with our mounted pumps. This semi-flexible bumpers help to isolate noise and eliminate creeping. For more rigid requirements, the pump base may be bolted directly to a firm foundation with or without the bumpers.

The pump should be located preferably in a clean and well ventilated area and adequate space should be provided whenever possible for routine maintenance such as changes of oil and belt adjustment and replacement. Above all, the pump should be located as closely as possible to its system in order to utilize it more efficiently.

Before connecting the pump, review the power source and the motor rating to be sure they agree in voltage, phase and frequency. On three phase applications the direction of rotation of the motor must be considered.

TYPES OF LUBRICANTS.

In order to obtain the maximum efficiency of your pump, it is important the use of Felisa Duo Seal oil (FE-1407K). Duo Seal oil has been especially prepared and is ideally suited for use in mechanical vacuum pumps because of its desirable viscosity, low vapor pressure and chemical stability. All Felisa vacuum pumps are normally tested with duo seal oil and shipped with a extra full charge of oil.

Duo Seal oil can be purchased with any of our dealers or directly with Felisa. The vacuum guarantee in our pumps applies only when Duo Seal oil is used.

TRAPS.

Where corrosive vapors or large quantities of condensable vapors are involved from vacuum processing, a trap may be used in the connecting line to the pump. It will help prevent damage to the pump mechanism and reduce oil contamination.

OPERATION.

a) Starting:

Before attaching the pump to a system it is well to familiarize yourself with the function and action of the pump which you have now acquired. Remove the intake and exhaust port plugs and temporarily provide a stopper for the intake and install the exhaust filter.

CAUTION: Do not run this pump with the intake open directly to the atmosphere. The pump will overheat, excessive oil mist will be emitted from the exhaust and the pump will eventually seize.

B) Oil level:

The amount of oil suitable for efficient and satisfactory performance should be determined after the pump has reached its operating temperature. If after a short period of operation the level should fall, it is likely the result of oil entering some of the interior pockets of the pump. Fill the pump until the oil level falls half way of the oil level window. If a gurgling sound occurs, additional oil must be added. Mechanical pumps will gurgle under three conditions of performance.

- * When the oil level in the pump is lower than required.
- * When a large leak is present in the system.
- * When the gas ballast is open.

Best performance of a mechanical pump is generally obtained after sufficient time has been allowed for the pump to come to operating temperature.

CLEANLINESS.

Take every precaution to prevent foreign particles for entering the pump. A fine mesh screen is provided for this purpose in the intake passage of all Felisa pumps.

LEAK DETECTION.

The importance of eliminating all leaks in a vacuum system is obvious when it is realized that a leak into the system, at atmospheric pressure, expands in volume by a factor of 750,000 to 10,000,000 or more. The pump must remove this added volume to maintain the desired vacuum. Fortunately a number of effective techniques for leak detection have been developed. Large leaks can be located by pressurizing the system and painting the suspected area with a tick soap solution. Escaping air will produce soap bubbles.

PRECAUTIONS.

Never use the pump for pumping fuels or fuel vapors an explosion might occur.

A few simple precautions are all that is necessary when shutdown is in order. If a gauge is connected to the system, first isolate the gauge, then turn off the pump and open the system to atmosphere. If the pump is removed from the system, cover the intake port with a rubber stopper or suitable cover to protect the pump against contamination and loose particles. If the pump has been contaminated in service and is going to be shelved for a prolonged period it is best to drain the oil and refill with a fresh charge.

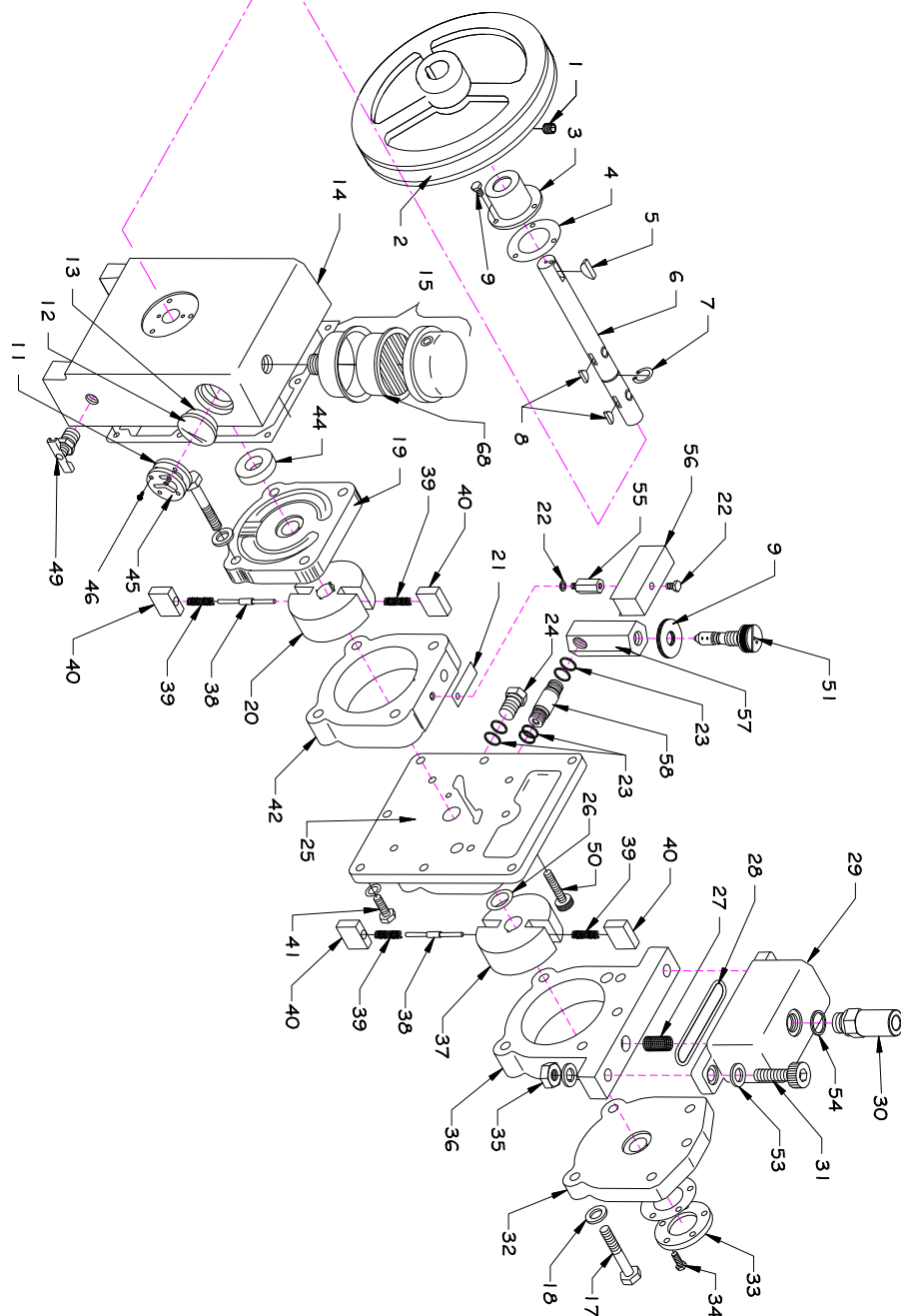
REPAIRS.

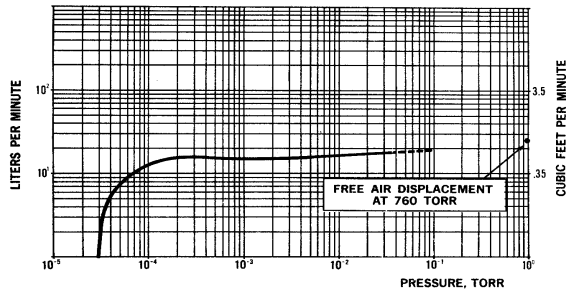
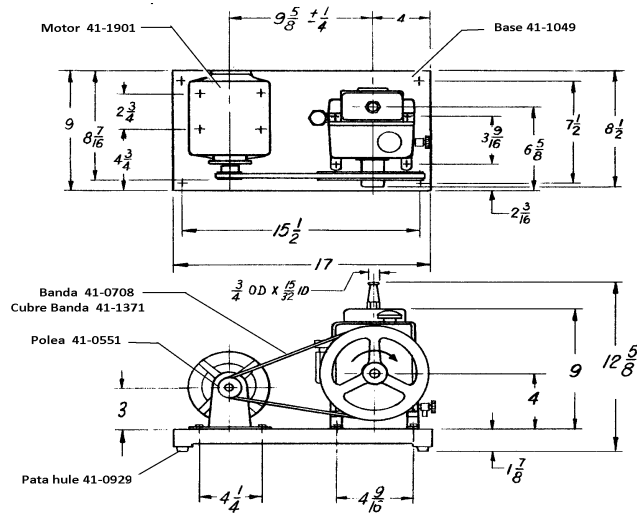
Like any other manufactured product, some parts of the pump could be damaged after a long operation time. To replace them always use original factory parts, a list of which is included in this manual. Always order the parts with its corresponding number. All the parts can be ordered from any distributor or directly from Felisa.

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You already have a one year warranty for your equipment. If after using this, answer this questionnaire and send it to us within 90 days by fax or email will automatically give one more year warranty!

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